



Anesthesia Shields

PROTECTS PATIENT

The shield provides a protective area around the patient's head and shoulders, reducing undue application of pressure to the patient and to monitoring lines.



DEFINES ANESTHESIA WORKING AREA

The shield provides a fixed and open space for the lines and vent tubes.



MOVES WITH THE OR TABLE

The shield easily mounts at the head of the OR table. Its feet slip into rectangular holders which mount to the OR table railing, allowing the shield to move with the table.





EXPANDS SURGICAL FIELD

The Pemco Anesthesia Shield has a tabletop surface that can be used for placement of surgical instruments, keeping them in the surgical field.



Anesthesia Shields

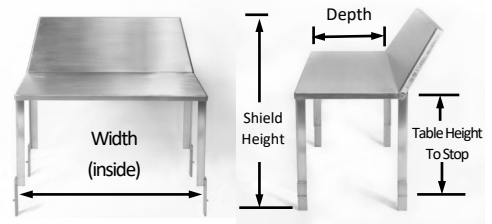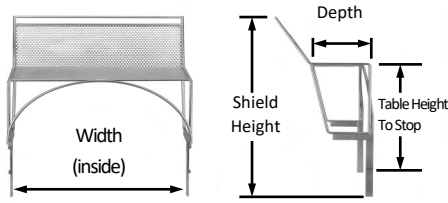

#4010-AC
Access 2

#4010-02
Dual Leg Anesthesia Screen

#4010
Four Leg Anesthesia Screen

Specs

Depth.....8.5 inches
 Width.....30 inches folded
 60 inches expanded
 Shield Height.....30 inches
 Weight.....12 pounds

Depth.....8 inches
 Width.....21½ inches
 Table Height..17 inches
 Shield Height.30 inches
 Weight.....12 pounds

Depth.....11 inches
 Width.....21½ inches
 Table Height..17 inches
 Shield Height.28 inches
 Weight.....,10 pounds

Features

Lightweight, strong, stainless steel tubing

Lightweight

Lightweight

Adjustable face shield

Slim profile, small shelf

Wide profile, wider shelf

Sturdy spring clips to hold drapes

Standardized flat posts for use with many OR tables.

Benefits

Protects patients facial area

Protects patients facial area

Protects patients facial area

Face shield is easily repositioned

Firm and sterile surface for surgeon to place instruments during surgery

Firm and sterile surface for surgeon to place instruments during surgery

Firm and sterile surface for surgeon to place instruments during surgery

Moves with the table

Moves with the table

Moves with the table

Tubing clips on the face shield allow the ventilator tubing and TEE probe to be supported, decreasing risk of injury to the patient

Drapes stay securely in place

Access II collapses for storage

Note the legs of the Dual and Four-Leg shield mount into rectangular bed rail holders and have permanent stops to prevent slippage.

Most OR tables are similar in width. Please confirm dimension of table and shield width. These must be within 1/2" in order for the shield to fit properly. Pemco can manufacture the shield to the required width if needed.

Contact Pemco or your Pemco representative for a quote.